

ENDOSCOPE INSTALLATION CHECK LIST

Dealership			
Contact		Position	
Telephone		E:	

Client			
Contact		Position	
Telephone		E:	
Address			

Building Type	Pre 1970's	Modern	New Build
Establishment Type	Private Surgery	NHS Surgery	Hospital
	Mobile Surgery	Training Establishment	Military
Ownership	Private	Franchise	Government
	Other	>	

No Surgeries at location	No Endoscopes required	
--------------------------	------------------------	--

CURRENT DENTAL SYSTEM

Dental Unit Brand			
Dental Unit Model			
Date Installed		Existing Warranty	Mnths
Dental Unit RAL Colour	<div style="display: flex; align-items: center; justify-content: center;"> <div style="width: 20px; height: 20px; background-color: #f08080; border: 1px solid black; margin-right: 5px;"></div> = Information for Fao </div>		
Internal Diametre of Unit Column			
External Diametre of Unit Column			
Height of Unit Column			
Existing Lamp at end of Column			

Compressed Air available?	Yes	No	To be installed
Power Available?	Yes	No	To be installed

Is a Monitor Already In Place?	Dental Class Monitor	15" Screen	17" Screen
	VGA Input Available	19" Screen	21" Screen
	Video Input Available	DVI Input	
Is a New Monitor Required?	Yes	No	Size (inches diagonal)
Is PC data required to be shown on monitor?	Yes	No	
Is there a requirement to capture images to PC?	Yes	No	Stills Video

ENDOSCOPE INSTALLATION CHECK LIST

ENDOSCOPE REQUIREMENTS

Variants in the INSIDE system are due to the length of the articulated arm and the type of installation (see notes - page 3)

	<input checked="" type="checkbox"/>	No.	Comments
Integrated Dental Unit Version - 750mm arm	<input type="checkbox"/>	<input type="checkbox"/>	<div style="border: 1px solid black; width: 100%; height: 100%;"></div>
Integrated Dental Unit Version - 400mm arm	<input type="checkbox"/>	<input type="checkbox"/>	
External Dental Unit Version - 750mm arm	<input type="checkbox"/>	<input type="checkbox"/>	
External Dental Unit Version - 400mm arm	<input type="checkbox"/>	<input type="checkbox"/>	
Ceiling Version - 750mm arm	<input type="checkbox"/>	<input type="checkbox"/>	

Integrated Applications	<input checked="" type="checkbox"/>	No.	
a. End Column Clamp Version	<input type="checkbox"/>	<input type="checkbox"/>	Code: 410600 (<i>Column diameter 45 or 50 mm</i>)
b. Clamp Version	<input type="checkbox"/>	<input type="checkbox"/>	Code: 410605 (<i>Column diameter 45 or 50 mm</i>)
c. Standard Version	<input type="checkbox"/>	<input type="checkbox"/>	Code: 410300 (<i>Column diameter 50 mm</i>)
d. Standard Version	<input type="checkbox"/>	<input type="checkbox"/>	Code: 410305 (<i>Column diameter 45 or 50 mm</i>)

External Applications	<input type="checkbox"/>	<input type="checkbox"/>	
a. End Column Clamp version	<input type="checkbox"/>	<input type="checkbox"/>	Code: 410660 (<i>Column diameter 45 or 50 mm</i>)
b. Clamp Version	<input type="checkbox"/>	<input type="checkbox"/>	Code: 410665 (<i>Column diameter 45 or 50 mm</i>)

Application Configuration	<input type="checkbox"/>	See Diagrams - page 3
---------------------------	--------------------------	-----------------------

Additional Notes / Comments

Installation Notes:





The INSIDE system can be fitted in two different ways on the dental unit:

- 1: Integrated Application: The main unit support needs to be changed.
Cables / connectors are inserted inside the lamp column
- 2: External Application: NO alteration Required.
Cables / connectors are outside the lamp column. In this case the cables are run inside a special protective conduit specifically created by Faro (see images - page 3)


- * Power Supply: 24v AC
- * Consumption: 0.7 A
- * Max Air Pressure: 6 Bar
- * Air Feeding Pipe: Internal Diameter 1.6mm / External Diameter 2.8mm
- * Air Consumption: 4 to 15 l/min

ENDOSCOPE INSTALLATION CHECK LIST

APPLICATION EXAMPLES

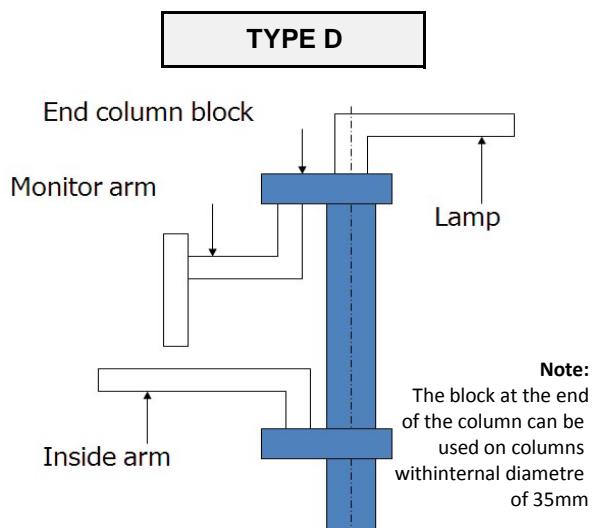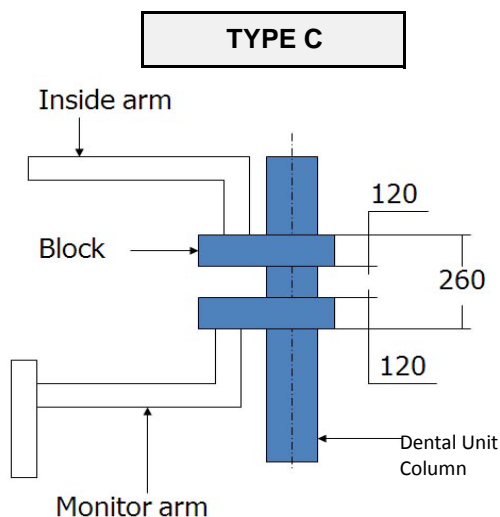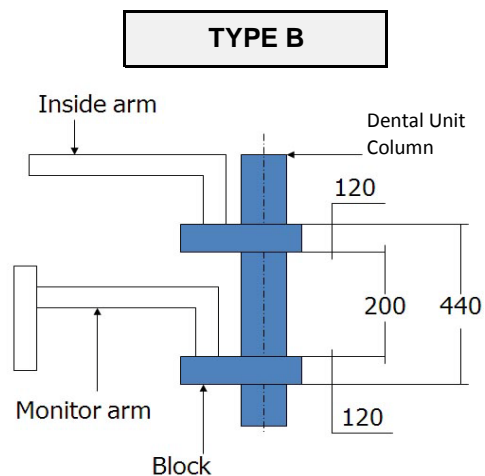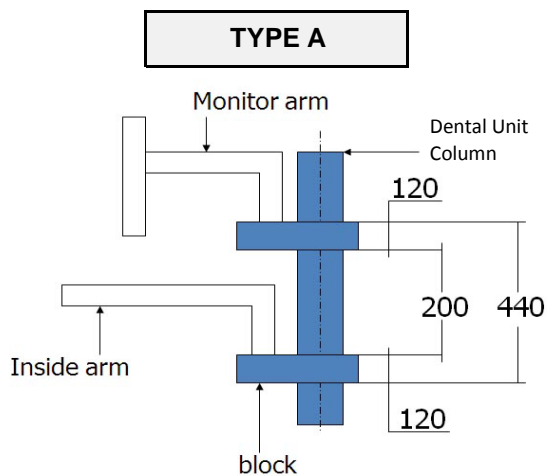

EXTERNAL APPLICATIONS

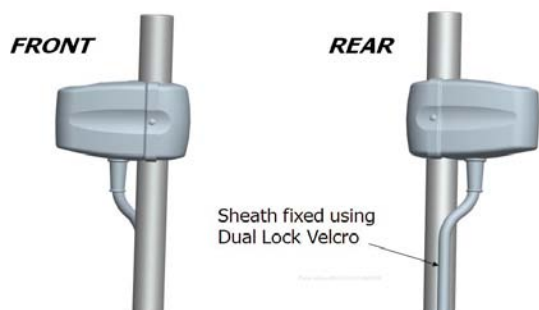

COLUMN CLAMPS

